


Securos Surgical Coxofemoral Luxation Management System

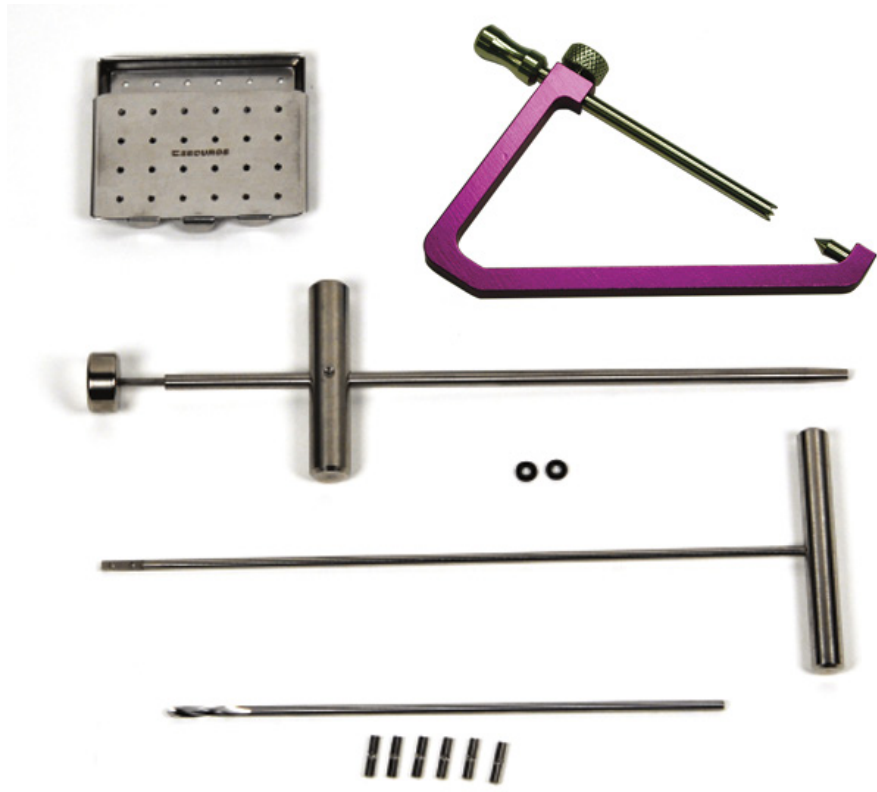

- Excellent alternative to a FHO for luxated hips.
- Toggle rods (pins) are made of the implant grade stainless steel.
- Tested and shown to be mechanically unrivalled.
- Aiming guide simplifies aiming and drilling of the femoral neck tunnel from the lateral side of the proximal femur.
- Insertion Device is used to place the toggle rod behind medial wall of the acetabulum.
- Suture puller is used to pull suture through drilled hole in femoral neck.
- The CLM System can be purchased as a kit or by independent components.

CLM Kit Includes:

- 1 - Femoral Head Aiming Guide
- 1 - Toggle Rod (Pin) Insertion Device
- 1 - Femoral Neck Suture Puller
- 2 - "O" Rings for Insertion Device
- 6 - Surgical Implant Toggle Rod (Pin)
- 1 - Small Sterilization Tray for Toggle Rods
- 1 - Orthopedic 3.2mm BITE Drill Bit

